

BIOHORIZONS[®]

single-stage



Single-stage Implant System

Science

Simplicity

Value

BIOHORIZONS[®]

BioHorizons® Single-stage Implant System

Specifically designed for quality outcomes

Simple Solutions™ restorations

Single Surgical Kit places BioHorizons Single-stage and Internal Implants.

1.8mm transmucosal machined collar helps establish biologic width.

Biomechanical square-thread design for high initial stability.

7mm implant length provides options in areas with limited bone height.

Both implant platforms and prosthetic components are color-coded for failsafe mating.


Simple Solutions Abutments can be restored without laboratory preparation.

Snap impression system requires no special tools.

PEEK Healing Cap for simple provisionalization.

Same proven hex connection as BioHorizons Internal Implant.

A full line of stock and custom castable abutments is also available.

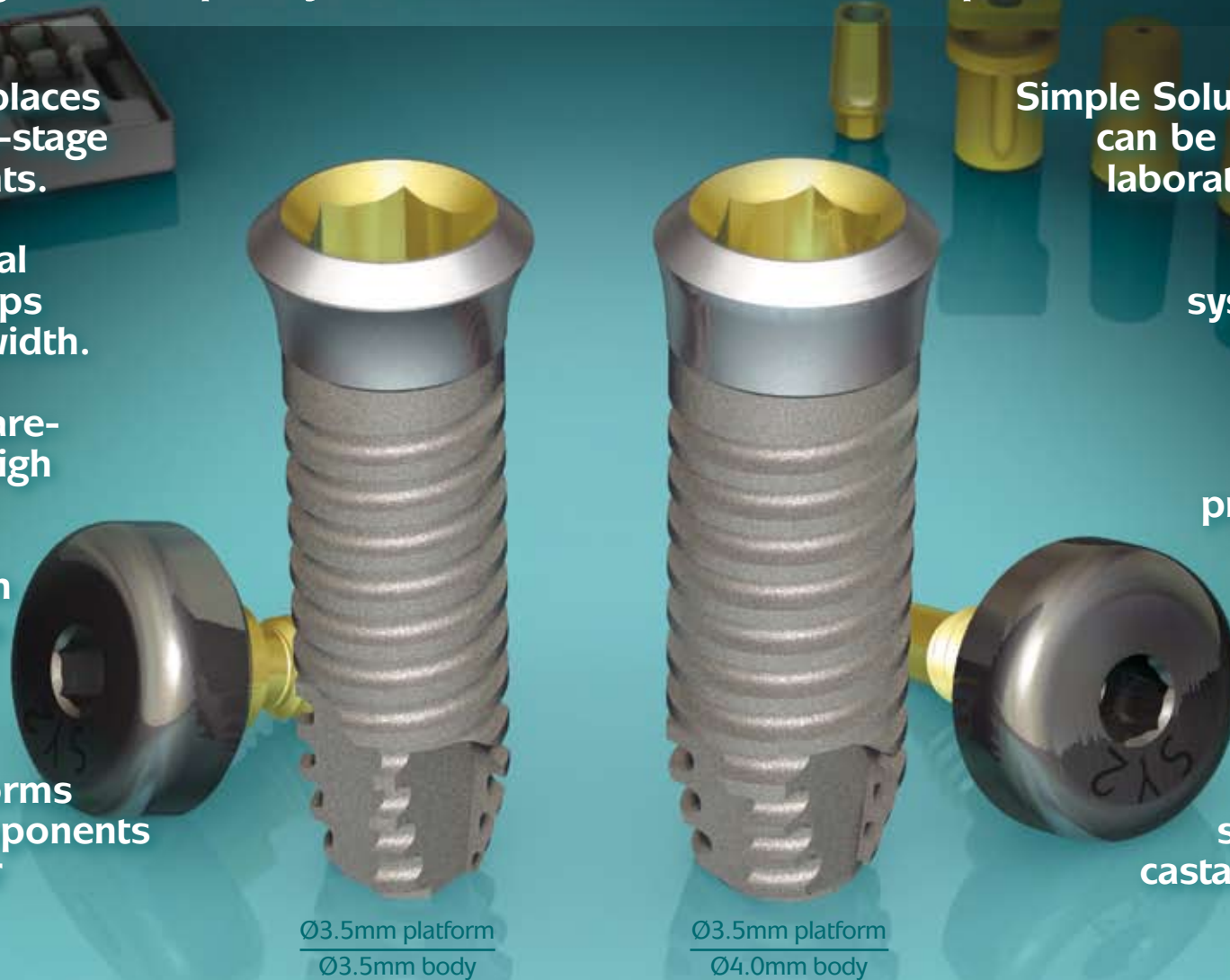

Ø3.5mm platform
Ø3.5mm body

Ø3.5mm platform
Ø4.0mm body


Simple Solutions Prosthetic Protocol

Each prosthetic platform offers two body diameter options

Seat Abutment

Make Impression

Deliver Prosthesis

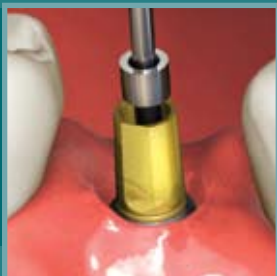

Prosthetic Platform

Ø3.5mm

Ø4.5mm

Ø5.7mm



Implant Body

Ø3.5mm | Ø4.0mm

Ø4.0mm | Ø5.0mm

Ø5.0mm | Ø6.0mm

A closed-tray pick-up of the Impression Cap captures the implant shoulder and Simple Solutions Abutment orientation. A provisional Healing Cap can be worn while the laboratory fabricates the final prosthesis.

Material: Titanium Alloy (Ti-6Al-4V)

Surfaces: Resorbable Blast Texturing (RBT) or Hydroxylapatite (HA)

Lengths: 7mm, 9mm, 10.5mm, 12mm & 15mm

BioHorizons • One Perimeter Park South • Birmingham • Alabama 35243

© BioHorizons Implant Systems Inc. 2007. All rights reserved.

**Call now for a complete informational packet on the
BioHorizons Single-stage Implant System**

or visit us on the web at

www.biohorizons.com

BioHorizons USA

Customer Care / Servicio al Cliente:

888-246-8338 or 205-967-7880

BioHorizons Canada

Customer Care / Service à la Clientèle:

866-468-8338 or 905-944-1700

BioHorizons UK

Customer Care:

+44 8700 620 550